

Polyclonal anti-human Erythropoietin

ORDERING INFORMATION

 Catalog No:
 hAP-0030

 Size:
 100 μg

 Storage:
 < -20° C</td>

Immunogen: human recombinant EPO

Ig G Type: Rabbit IgG

Applications: Western blot & Immunohistochemistry

Preparation: Produced in rabbits immunized with purified recombinant human erythropoietin (rhEpo).


Formulation and Storage: The anti-serum was purified by Protein A/G affinity chromatography and lyophilized from a 0.2 μm filtered solution in phosphate-buffered saline (PBS). Lyophilized samples are stable for 2 years from date of receipt when stored at -70°C.

Reconstitution: Reconstitute the antibody with 500 μ l sterile PBS and the final concentration is 200 μ g/ml. Reconstituted antibody can also be aliquoted and stored frozen at < -20 for at least for six months without detectable loss of activity. **Avoid repeated freeze-thaw cycles.**

Specificity: This antibody recognizes human EPO in western blot. Reactivity to EPO of other species is not tested.

Applications:

Western blotting: Yes (1:600-2000)
 Immunohistochemistry Yes (1:50-200)
 (Paraffin sections)

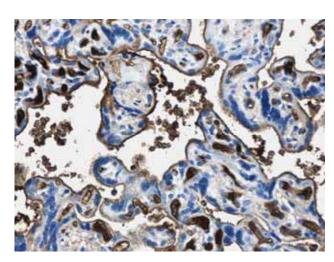

Formalin fixed and Paraffin embedded human placenta tissue section was subjected to IHC staining using hAP-0030.

Contact & Ordering Information: Angio-Proteomie, 11 Park Drive, Suite 12, Boston, MA 02215, USA. Tel: 001-6175492665; Fax: 001-4802474337, angioproteomie@gmail.com