

Mouse Monoclonal anti-human EDG-1

ORDERING INFORMATION

Catalog No:	hAP-0045
Size:	100 µg
Storage:	< -20° C
Immunogen:	human recombinant EDG-1 protein
Ig Type:	mouse IgG1
Clone:	MM0045-21L9
Applications:	FC, Western Blot and IHC

Preparation: This antibody was produced from a hybridoma (mouse myeloma fused with spleen cells from a mouse immunized with human recombinant Endothelial Differentiation Gene-1 (EDG-1) protein.


Formulation and Storage: The IgG fraction of culture supernatant was purified by Protein G affinity chromatography and lyophilized from a 0.2 µm filtered solution in phosphate-buffered saline (PBS). Lyophilized samples are stable for 2 years from date of receipt when stored at -70°C.

Reconstitution: Reconstitute the antibody with 500 µl sterile PBS and the final concentration is 200µg/ml. Reconstituted antibody can also be aliquoted and stored frozen at < -20 for at least for six months without detectable loss of activity. **Avoid repeated freeze-thaw cycles.**

Specificity: This antibody recognizes human EDG-1 in western blot. Cross-reactivity to other species has not been tested.

Applications:

1. Western blotting:	Yes	(1:100-1000)
2. Flow Cytometry	Yes	(1:50-200)
3. IHC (Formalin/Paraffin)	Yes	(1:50)

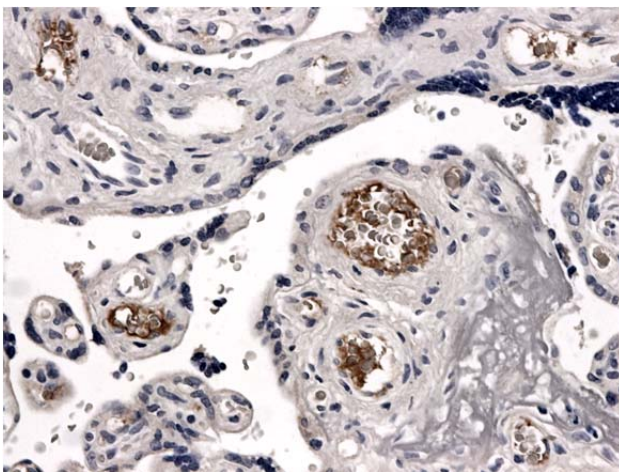

Formalin fixed and paraffin embedded normal human placenta tissue is subjected to mouse anti human EDG-1 (hAP-0045)

Final dilutions should be optimized by the end users and the dilutions provided in the data sheet are for reference only.

Contact & Ordering Information: Angio-Proteomie, 11 Park Drive, Suite 12, Boston, MA 02215, USA. Tel: 001-6175492665; Fax: 001-4802474337, angioproteomie@gmail.com