


## Product Specification

<b>Product No.</b>	HGPX2-05 (-1, -2.5, -5)
<b>Lot No.</b>	HGPX2210722
<b>Product Description:</b>	Recombinant Human Glutathione Peroxidase 2 (Non-tagged; 190 residues)  0.5 (1, 2.5, 5) Units, ca 50 (100, 250, 500) µg protein stored in 25 mM Tris-HCl (pH 8.0) with 250 mM KCl, 2 mM 2-Mercaptoethanol and 10% glycerol
<b>Specifications:</b>	1 U of enzyme consumes 1 µmol GSH per min in a Glutathione Reductase-coupled GPx assay containing 1 mM GSH, 0.5 mM Cumene Hydroperoxide, 15 nM Human Glutathione Reductase and 0.3 mM NADPH in TE buffer (50 mM Tris-HCl, 2 mM EDTA, pH 7.5)
<b>References:</b>	Cheng Q, Arnér ES. (2017) Selenocysteine Insertion at a Predefined UAG Codon in a Release Factor 1 (RF1)-depleted Escherichia coli Host Strain Bypasses Species Barriers in Recombinant Selenoprotein Translation. J Biol Chem. 292:5476-5487.
<b>Storage:</b>	Shipped in cold pack by courier service. Store at -20° C upon delivery.
<b>Examples of Use:</b>	Glutathione Reductase-coupled GPx assay; Inhibitor screening; etc.
<b>NOTE!</b>	Centrifuge the tube briefly before opening the lid
<b>Purity and Gel Analysis:</b>	Purity is ≥ 95% as judged by gel analysis

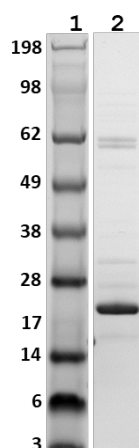

*Coomassie stained SDS PAGE analysis of 30 µg HGPX2-50 product*  
*Lane 1: Molecular weight marker (SeeBlue™ Plus2)*  
*Lane 2: HGPX2-50*